

Never hassle with shrink bands again!



The Weber DB™ Sterilized Pre-filled Dilution Bottle has a sterility strip attached to the cap. Much easier to open than shrink-wrap, especially when wearing gloves.

DBBR1111

Copyright © 2011 Weber Scientific

WEBER SCIENTIFIC

2732 KUSER ROAD, HAMILTON, NJ 08691

Please let us know

- ✓ any corrections or changes to your address
 - ✓ if this name should be deleted from our list
- We need your help – thanks!

800-328-8378

www.weberscientific.com

PRSRT STD
U.S. POSTAGE
PAID
U.S. Press
31601

DB™ Pre-filled Dilution Bottles



Don't hassle with shrink bands!




Easy-pull
sterility strip



Easier to store
Greener packaging
60% smaller footprint and 3 times less
corrugated waste than the competition

Much easier to open, use and store!

WEBER SCIENTIFIC

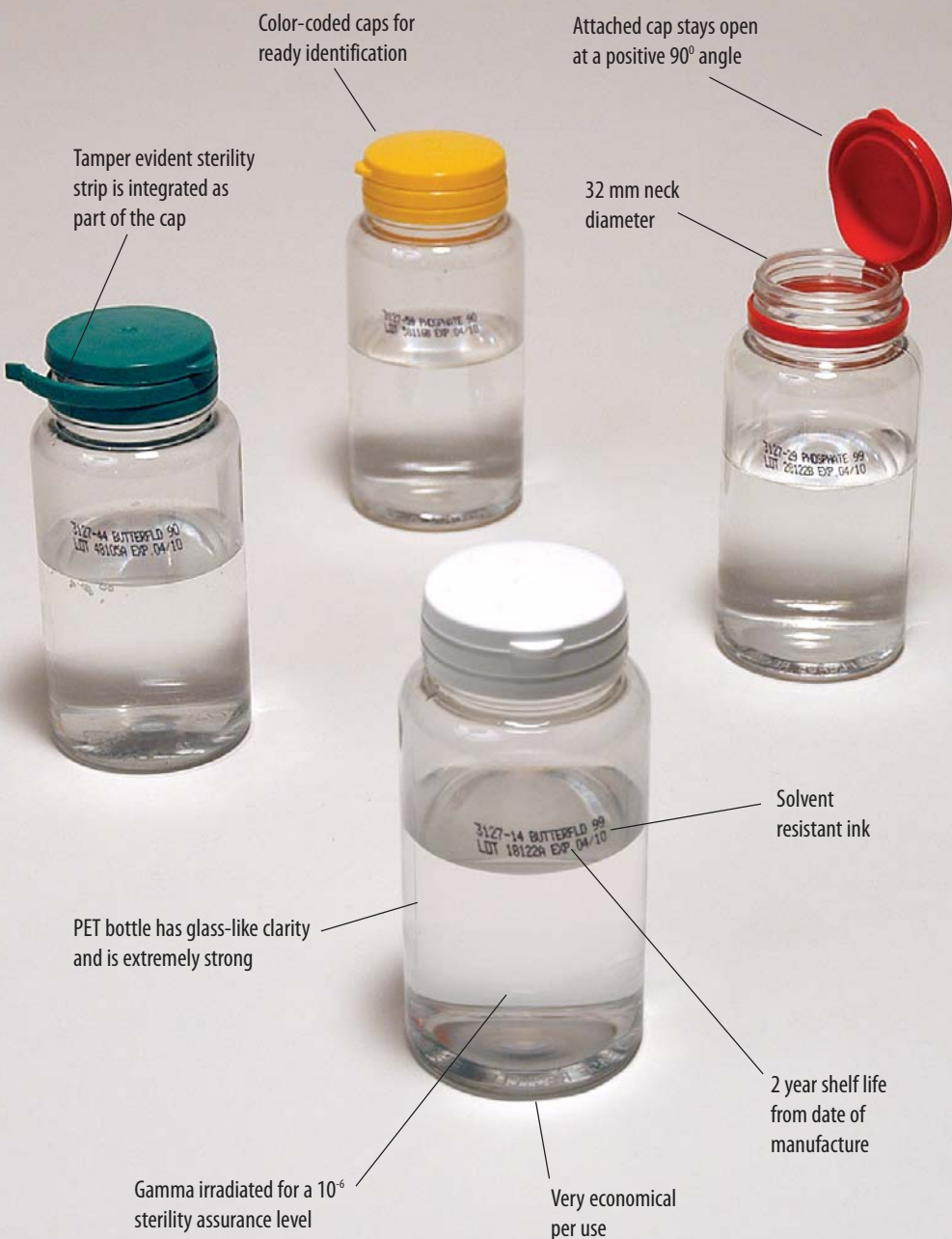

Butterfield's Buffer for Foods and Dairy Products (APHA, FDA, AOAC or USP methods)
Contains purified water, monobasic potassium phosphate, and sodium hydroxide for pH adjustment.

Phosphate Buffer for Water/Wastewater and Dairy Products (APHA or EPA methods)
Contains magnesium chloride, monobasic potassium phosphate, and sodium hydroxide for pH adjustment.

Every bottle is pre-filled to either 99 or 90 mL ± 2 mL, buffered to pH 7.2 ± 0.2, has an overall capacity of 150 mL, and is recyclable.

Catalog Number	Formulation	Fill, mL	Cap Color	Available	Your Cost per case	Your Cost each
3127-12	Butterfield's	99	white	} case/72	1-11 42.84	59.5¢
3127-44	Butterfield's	90	green			
3127-29	Phosphate Buffer	99	red			
3127-59	Phosphate Buffer	90	gold			
					12-24 37.44	52¢
					25+ 35.28	49¢

Ordering Hint: Our case contains 72 bottles, compared to 66 for some brands.

Large Volume Pricing is Available!

Do you use a lot of dilution bottles? We'd be delighted to quote high volume pricing, over and above our heavily discounted 25 case price break. Please contact your Account Manager.



Sharon Wilson
National Accounts Manager
swilson@weberscientific.com



Nancy Silvester
Account Manager
nsilvester@weberscientific.com



Danielle Spahn
Account Manager
dmspahn@weberscientific.com

FEATURE / BENEFITS

ADVANTAGES

Lot # and exp. date printed on side of bottle

solvent resistant ink

Confidently sterilize outside of cap

- ink withstands 70% Isopropyl Alcohol

PET bottle

glass-like clarity
exceptional strength

The only PET disposable dilution bottle available

- permits meticulous visual inspection of contents
- virtually eliminates cracking or breaking in handling or shipping
- safer - avoid release of contents if dropped

Gamma irradiated

10⁻⁶ sterility assurance level

The highest sterility level available

- no chance of microbial contamination - none

Attached cap

no cap to set down or hold

- enhances aseptic handling and reduces chance of contamination
- promotes use of two hands for superior ease-of-use
- only one sample identification number needed, not two
- stays open at a positive 90° angle
- leak-proof - samples can be vigorously shaken without fear

Sterility strip attached to cap

tamper evident

You'll love not having to hassle with shrink-bands

cap stays securely sealed

- much easier to open than shrink-wrap, especially when wearing gloves
- can be shipped by ground or air - no transportation limits - holds up to rough handling

Color-coded

rapid and accurate identification

- recommended when using multiple formulations or fills
- advances good laboratory practice

32 mm neck diameter

wide-mouth

- convenient for liquid or solid samples

2-year shelf life from date of manufacture

extended expiration dates

The longest shelf-life in the industry

- no worry about expired product; no waste saves money

This product is made in strict compliance with specifications as published in the most recent edition of A.P.H.A. Standard Methods for the Examination of Dairy Products; A.P.H.A. Standard Methods for the Examination of Water and Wastewater; A.P.H.A. Compendium of Methods for the Microbiological Examination of Foods; F.D.A. Bacteriological Analytical Manual; F.D.A. Milk Laboratory Evaluation Form 2400; United States Pharmacopoeia.

Each lot is thoroughly tested in house with toxicity done by an independent testing laboratory, and is released only after passing all specifications for pH, sterility, toxicity, fill volume, and package integrity and markings. The COA indicates that each particular lot passes all of these specifications. Specific results are given for pH and toxicity analysis.

You may obtain a copy of the Certificate of Analysis (COA), by lot number, online at weberscientific.com.

TEST	RESULTS	SPECIFICATION	METHOD
pH @25°C	7.29	7.2 -- 8.2	WS-00P-8.2
TOXICITY*	Clear at 0 minutes; Clear at 30 minutes; PASSES	* 20% DECREASE	Perkin Elmer 1122M-0906
STERILITY	PASSES		Perkin Elmer 1122M-0906
FILL VOLUME	CONFORMS	± 0.10 ml	ASAO - VALIDATED (Commercial)
COLOR	CONFORMS	LABEL ± 2 MI	WS-00P-8.2, 11.2
PACKAGE INTEGRITY STRIP*	CONFORMS	CLEAR NO FOREIGN MATTER	VISUAL WS-00P-11.2
PACKAGE INTEGRITY BOTTLE	CONFORMS	NO VISIBLE DAMAGE TO CARBON	VISUAL WS-00P-11.1
PACKAGE MARKINGS	CONFORMS	NO LEAKING	WS-00P-14.1
BOTTLE MARKINGS	CONFORMS	CONFORMS TO STANDARD	WS-00P-8.1

*TEST PERFORMED AT COMMERCIAL LAB. (Perkin Elmer, Inc.)
**Approved by the Association of Official Agricultural Chemists (AOAC) method 4500.01-0101
***ISO 11731:2003 (ISO 11731:2003) method AAM (2003)

Made in:
Weber Scientific
1122 River Road
Harrison, NJ 07033
(908) 584-7670

Paul Weber
President

Date: 08/06/14 (2014)

Would you like to evaluate a free sample?

No problem . . . let us know and free samples will be sent right out! Please advise which of the four products you would like to evaluate when making your request.



ORDERING INFORMATION

How to Order



PHONE:
800-328-8378

8 a.m. to 5:30 p.m. EST
(or 609-584-7677)



FAX:
609-584-8388



MAIL:
WEBER SCIENTIFIC

2732 Kuser Road
Hamilton, NJ 08691-9430



ONLINE:
weberscientific.com

E-MAIL:
info@weberscientific.com

✓ **No Minimum Order.**

✓ **In-Stock for Immediate Shipment.
We Mean It!**

✓ **Order Risk-Free with Our No-Nonsense
Satisfaction Guarantee.**



If you are not totally satisfied, call our customer service department within 30 days for replacement or return authorization.

Huge freight savings and faster delivery throughout the United States

Dilution bottles now ship out of our Nevada warehouse for delivery to Western states and out of our Missouri warehouse for selected Southern and Midwestern states. Customers in these areas will realize significant freight savings of at least \$1.00 to over \$10.00 per case. Average volume shipments will save over \$6.00 per case.

Plus, Weber DB™ shipments will be received between two and five days sooner.

Eastern states will continue to receive fast and economical service from our facility in New Jersey.

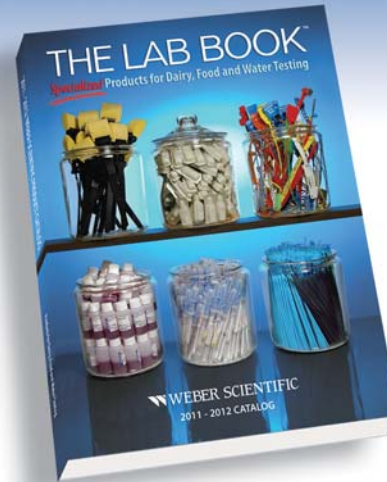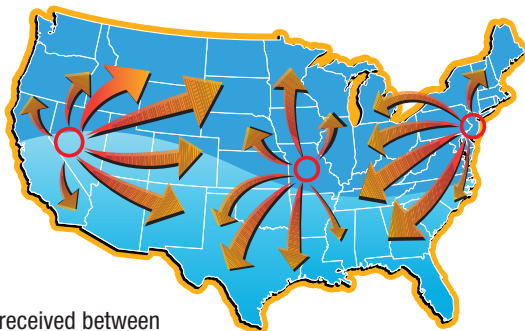

FREE

Weber Scientific NEW 2011-12 Catalog

Be certain to request your FREE copy of this essential 265-page Buyers Guide dedicated to Dairy, Food and Water Analysis Equipment and Supplies.

**Uniquely arranged by Quality Control
Test category ...**

- Apparatus
- Bacteriological Testing
- Books
- Chemicals
- Dairy Testing
- Juice Testing
- Labware
- Liquid Handling
- pH & ION Determination
- Refractometers
- Safety Supplies
- Sampling Supplies
- Sanitation Supplies
- Security
- Thermometers
- Timers
- Water and Wastewater Testing

**It's easy to request
your FREE copy ...**

Call **800-328-8378** 8 a.m. to 5:30 p.m.
Eastern Time (or 609-584-7677)

Or go online at **www.weberscientific.com**



- ✓ **Weber Scientific is a laboratory supplier that you will actually enjoy doing business with.** More than ever we hear that the mega-sized suppliers aren't particularly focused on the specialized needs of food and dairy labs, and that their level of service is often indifferent.
- ✓ **Great Product Selection.** We're specialists in quality control supplies for food and dairy labs. You'll find hundreds of laboratory essentials plus many exclusive and hard-to-find items.
- ✓ **Legendary for great prices.** We offer some of the best prices in the industry without any of the dickering or hassles and without committing to any complicated or long-term contracts.
- ✓ **Free freight upgrades.** Receive expedited service nationwide at our expense – not yours.
- ✓ **100% Satisfaction Guarantee.** Count on it! We've worked hard for over 50 years to earn your trust.
- ✓ **Telling our Customers Thank-You – again and again.** Our *Frequent Buyers* Program is a way for us to express our continued gratitude for your repeated support.

Order Toll-Free **800-328-8378**

Visit **weberscientific.com**

 **WEBER SCIENTIFIC**